meodental

Meodental Ltd.

Breakthrough Technology in Teeth Whitening
& More

meodental

While the process of decomposition of H_2O_2 depends on :

- Temperature
- Impurities
- Stabilizers

And also :

- Metals of transition
- Catalyst

The **Natural+™** patented, ecological active agent has a literally different operating mode.


 **NATURAL+™**
The natural solution for teeth whitening

meodental

Since the **Natural+™** active agent is stabilized, its process of releasing oxygen is *progressive* and not volatile and its whitening power is increased.

This specific property allows the formula to keep its efficiency over a period of more than **30 min.** thus *avoiding the need to renew the application of the gel between two sessions of treatment.*

Whitening efficiency

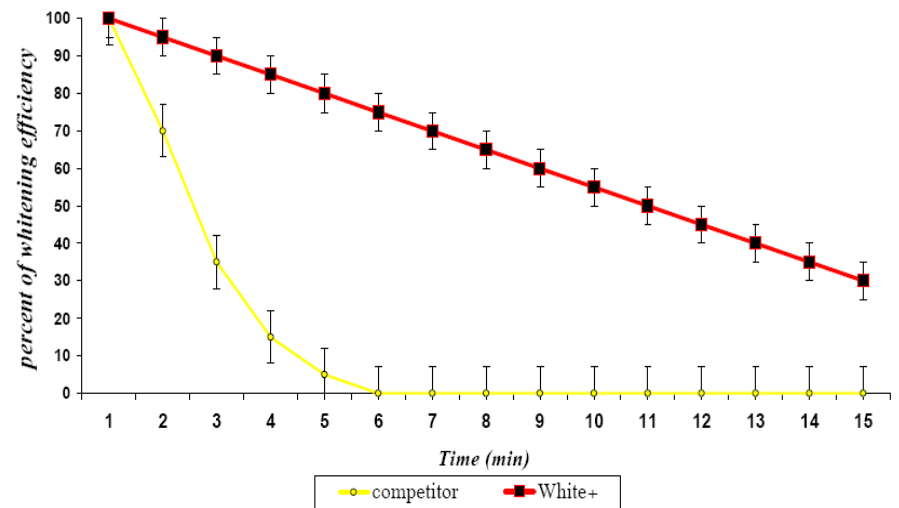NATURAL+™

The natural solution for teeth whitening

meodental

The light of the **Natural+™** 4th generation LED Lamp, without giving off heat or UV radiation, controls the whitening properties of the **Natural+™** Gel.

The **Natural+™** LED Lamp is a revolutionary LED (Light Emitting Diode) system and is the only one on the market that has several scientifically developed programs designed to accelerate and optimize the results of the **Natural+™** Gel as well as any other dental whitening product.

The blue LEDs act as a catalyst to the whitening gel and improve the efficiency of whitening superficial and in-depth tooth discolorations.

The red LEDs assist in reducing sensitivity when working with the **Natural+™** Gel unique formula.

Thanks to the stability of the formula, the dentist can control the pace of decomposition via the intensity and duty cycle of illumination by selecting the appropriate **Natural+™** LED Lamp program.

And now, with the higher intensity, the **Natural+™** LED Lamp can also operate as a curing light for polymerizing the gingival protection barrier.

 **NATURAL+™**
The natural solution for teeth whitening

meodental

The **Natural+**TM Gel includes *Calcium* in the formula:

- To prevent against demineralization
- To act as a buffer on the acidity
- To allow the whitening gel to preserve a *neutral pH*

NATURAL⁺TM

The natural solution for teeth whitening

meodental

Since with **Natural+™** you reach a greater whitening power, it significantly reduces the treatment time at the dentist's chair.

A single session of 30 min. is enough to gain up to 8 shades.

NATURAL+™

The natural solution for teeth whitening

The **Natural+**TM gel is enriched in:

- *Potassium Nitrate*, which role is to affect the threshold of reaction of the sensory nerves and so to prevent pain
- *Sodium Fluoride*, which strengthens the calcic barrier enamel / dentin and has a cariostatic action

Healthy Smile or Hidden Damage?

Images adjusted to illustrate clinical situation

meodental

The **Natural+**™ formula maintains an optimal degree of hydration of the tooth all through the dental bleaching procedure, thus avoiding the exchanges of oral flows.

NATURAL+™

The natural solution for teeth whitening

meodental

Meodental's **Natural+™** is a unique process which answers even the most difficult clinical cases as the treatment of *dyschromies of intrinsic type (tetracycline / fluorosis) and extrinsic type*.

 **NATURAL+™**
The natural solution for teeth whitening

Extrinsic dyschromias

Before
A3/A4

After
B1/B0

Intrinsic dyschromias

Before
A3/A4

After
B0-SB3

Decalcification

Before
A3/A4

After
B1/B0

Grey tooth & Tetracycline stains

Before
A3/C3

After
A1/C1

Fluorosis

Before


After

meodental

Natural+™ Product Range

The **Natural+™** Professional Kit is intended for in-clinic application by a Dentist.

It contains our exclusive **Natural+™** formula, our highly efficient and safe to use whitening solution plus all the required ingredients to perform the complete procedure and obtain effective and lasting results.

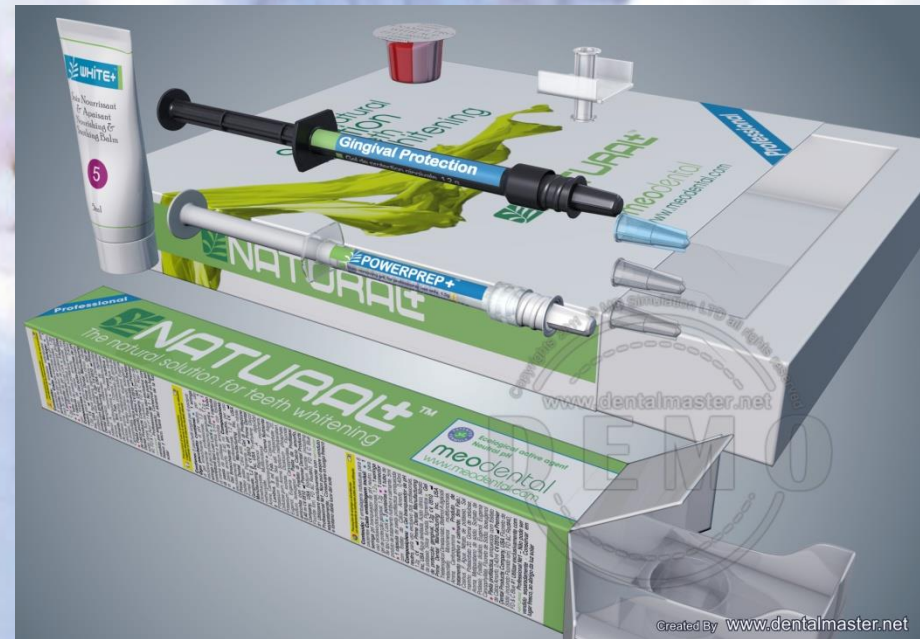 **NATURAL+™**
The natural solution for teeth whitening

meodental

Natural+™ Product Range

The **PowerPrep+™** Professional preparation and cleaning gel intended for in-clinic application by a Dentist prior to the bleaching process in order to reduce treatment duration, improve results and reduce sensitivity..

The PowerPrep is a unique formula combining citric acid for:

- Cleaning superficial stains.
- Removing smear layer.
- Opening tubules (speeds bleaching).
- Assists in dissolving de-calcification stains.

And potassium nitrate for:

- De-sensitizing the sensory nerves.

The main advantage of the combination is that it has a neutral PH and, therefore, does not etch the teeth.

NATURAL+™

The natural solution for teeth whitening

meodental

Natural+™ Product Range

Natural+™ Home Kit is a self administering kit that patients can obtain from the dental clinic after consultation.

It ensures optimal results from the comfort of the patient's home as a standalone treatment or, for more severe cases, in combination with the **Natural+™** Professional application.

It also contains our exclusive **Natural+™** Gel which is safe, reliable, and far more efficient than any other home treatment or over the counter product, without the side effects or sensitivities.

Meodental's unique constant-dose application syringe delivers the correct amount of gel per tooth thus ensuring patient safety , comfort and ease of use.

 **NATURAL+™**
The natural solution for teeth whitening

meodental

Natural+™ Product Range

Meodental's **PowerPrep+™** toothpaste is used in preparation before home whitening procedures.

Its formulation allows the cleaning of unwanted on-surface stains.

The presence of desensitizing agents in the formula prevents against dental sensitivities.

It:

- Prepares the teeth without affecting them
- Optimizes the results obtained by the teeth bleaching treatment
- Speeds the whitening process
- Minimizes sensitivities
- Its pH preserves the oral balance.

 **NATURAL+™**
The natural solution for teeth whitening

Differences (1/2) – Whitening gel

- First and foremost the hydrogen peroxide used in our gel is totally different, mainly because it is stabilized (patented).
If regular hydrogen peroxide dissolves in an open bottle at room temperature within 3-5 minutes ours is stable for up to 3 YEARS.
- In addition, thanks to the different manufacturing process of our gel, it does not include any residues of heavy metals.
These residues, present in all other products, connect to minerals on the tooth surface and cause diffusion of water from the tooth and, thus, dehydration.
- Another advantage of our stabilized formula is in our ability to add a high concentration of minerals without dissolving the hydrogen peroxide.
- Our formula, unlike others, maintains a neutral PH not only in the syringe but also on the tooth while activated and throughout the treatment.

Natural+™ – Differences & Advantages – Summary

Differences (2/2) - Supporting products

- The patented formula of the **PowerPrep+**™ in both its professional gel form or home toothpaste form prepare the teeth for the whitening while assisting in reduction of sensitivities.
- The #5 re-mineralizing and soothing balm, supplied in both the professional and home kits, resolves residual post-treatment sensitivities while strengthening the tooth by the addition of healthy minerals.
- The **Natural+**™ LED Lamp controls the activation of the gel without heating it through illumination at 470nm, and reduces sensitivity through the use of red light at 660nm. In addition the new and improved ergonomic design along with the audio ability increase patients' comfort and calmness during the procedure.

Advantages (1/3)

- Virtually no sensitivity –
 - The gradual linear (rather than fast exponential) release of the radical oxygen ions.
 - Neutral PH throughout the treatment.
 - No tooth dehydration
 - High concentration of minerals with desensitizing effect as calcium and potassium nitrate.
 - The use of **PowerPrep+™** prior to the treatment and #5 balm after it resolve even the rare cases of remaining sensitivities.
 - The red light in the lamp assists in reducing sensitivities particularly in the gingiva.

Advantages (2/3)

- Simpler procedure –
 - Only one application is required for the entire in-office treatment of 30-45 minutes with gel and lamp.
 - Entire procedure (chair-time) duration of less than one hour (about half of any other product).
 - In all regular cases no second session or home continuation is required.
 - Both arches are treated simultaneously without the risk of sensitivity.
 - Severe cases, previously considered untreatable, (e.g. tetracycline stains) are treated with the exactly same simple procedure only with either home continuation or a second session.
 - All the required “supporting” products are supplied in the kits simplifying the dentist’s preparation.

Natural+™ – Differences & Advantages – Summary Advantages (3/3)

- Higher efficiency –
 - **PowerPrep+™** use prior to the whitening gel removes superficial staining and speeds up whitening.
 - The intensity of the blue light in the lamp is optimized for whitening with the **Natural+™** gel
 - Continuous slow release of the radical oxygen ions enables a larger portion of them to penetrate the tooth through the enamel.
 - Tooth is not dehydrated (rather is constantly irrigated) so the oxygen ions, carried by water, penetrate deeper into the dentine.
- Long-lasting results –
 - No regression due to tooth re-hydration (there is no dehydration to begin with). Teeth maintain transparent and vivid look after the treatment.
 - Slow release of oxygen ions continues up to 24 hours after treatment ends resulting in further improvement.

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A progress in the practice of professional bleaching technical at the clinic

Clinical case1:

- Picture1: Photo after cleaning and before bleaching.
Note the defect of the enamel on the 22.
- Picture2: After treatment (with no pain).
Note the penetration of the product in the defect of 22.
- Picture3: Shade is stabilized after one week (cigarette smoker).
The bleaching effect is visible on the 22 which is taken here as a reference.

Clinical case2:

- Picture4: After cleaning and before whitening.
Starting shade 2M 1, Vita scale 3D Master.
- Picture5: The luminosity is apparently increased.
Checking with starting shade referenced 2M1.
- Picture6: New checking with a clearer shade, 1M 1 showing the whitening effect.
- Picture7: Checking after 24 hours with this technique, stabilized results.
- Synthesis: These two clinical cases reveal that the White+™ bleaching system from Meodental is reliable.
The ergonomic equipment is compact and the procedure is simplified.
There are less sensitivities (especially during the treatment) and the results are stable.

Dr Herszenfis Frederic: Master's degree in Biological and Medical Science, Partner of the SELARL Touati-Herszenfis.

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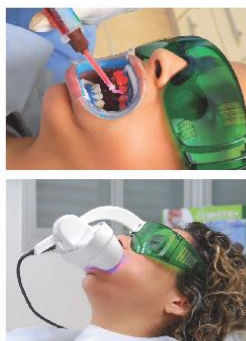
White+™: Revolutionizing teeth whitening

All the dental bleaching products on the market which contain hydrogen peroxide work on the same mode, because the process of decomposition of the hydrogen peroxide applies to all on the same basis.


In brief, its process of decomposition is completed in an exothermic reaction of dismutation into water and dioxygen, its rate being influenced by factors depending on the temperature, the impurities and the stabilizers. Commonly, it is metals of transition and derivatives which facilitate the decomposition of the hydrogen peroxide. The use of a catalyst also comes to increase strongly its decomposition into oxygen and gives off heat.

White+™ is a patented, ecological active agent, which prevents these phenomena. Its operating mode is literally different, and its whitening power is increased, safely.

It works on an "atypical" mode because, stabilized, its mode of decomposition is progressive and not volatile, as inherently hydrogen peroxide is.

It presents numerous advantages. It is not necessary to renew the pose of the gel between two sessions to reactivate the efficiency of the product which releases the oxygen over a period of more than 20 minutes.

Whitening efficiency

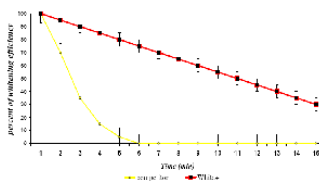It does not oxidize because it includes in its composition chelator agents which allow it to keep its stability. It therefore departs from all the common process, existing on the market, as it reacts according to its own operating mode. Thus it is not necessary anymore to bring it together with metals of transition (ex. Zoom2) in order to potentiate its activity. The light of the White+™ lamp, without giving off heat, allows activating its whitening properties, without producing residues. Since White+™ is stabilized, the dentist can control the pace of decomposition via the intensity and duty cycle of illumination by selecting the appropriate White+™ lamp program.

The active agent White+™ includes Calcium at an optimal rate to prevent against demineralization of the tooth and while the hydrogen peroxide is entering into its phase of activation, it reacts on the Calcium, included into its composition, which acts as a buffer on the acidity and then allows the bleaching gel to preserve a neutral pH.

Prolonged exposure of teeth to bleaching compositions, as practiced at present, has a number of adverse effects in addition to that of tooth sensitivity. These include: solubilization of Calcium from the enamel layer at a pH less than 5.5 with associated demineralization; penetration of the intact enamel and dentin by the bleaching agents, so as to reach the pulp chamber of a vital tooth thereby risking damage to pulpal tissue; and dilution of the bleaching compositions with saliva with resulting leaching from the dental tray and subsequent digestion.

The exceptional properties of the active stabilized agent White+™ and the presence of Calcium in the formula, allow when these two components are associated to release the maximum of oxygen in a constant and progressive way during the duration of the treatment. As a result it has a greater whitening power and a reduction of the processing time at the dentist's chair for the patient (a single session of 30 min. is enough to gain up to 8 shades).

The White+™ gel, besides containing Calcium is enriched in Potassium Nitrate, which role is to affect the threshold of reaction of the sensory nerves and so to relieve the pain as well as Sodium Fluoride, which strengthens the calcic barrier enamel/dentin and acts as a desensitizing substance of has a cariostatic action.

Tooth sensitivity is believed to result from the movement of fluid through the dentinal tubes toward nerve endings in the tooth. This movement is enhanced by the carriers for presence of components which can be found in other bleaching gel compositions.

Meodental's White+™ is a unique process, benefiting of a neutral pH, which answers to all clinical cases, even the most difficult as the treatment of dyschromies of intrinsic type, among which those bound to antibiotics (tetracycline) and fluorosis. White+™ is also very effective for the treatment of the age-related extrinsic dyschromies, or in case of consumption of tea, coffee, wine, cigarettes, etc.

The White+™ formula allows maintaining an optimal degree of hydration of the tooth during all the procedure of dental bleaching, so avoiding the exchanges of oral flows.

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Professional Technical at the clinic

- Technique : Zoom
(Discus Dental – Laboratory Pierre Roland)

Picture 3a: 37 years-old patient - Uniform yellow shade (B3 on Vita scale Classic). We used a combined technique. First, we realize a treatment at the dentist's chair, with the Zoom lamp (treatment duration= 1h20min) followed with six nights of home treatment (carbamide peroxide at 10% of concentration) to stabilize the results. Preparation and isolation of the teeth are essential to guarantee the success of this session.

Picture 3b: Results at the end of treatment (dentist's chair + home), one week of treatment all in all.

- Technique : White+™
(Meodental)

Picture 4a: Starting case: Twenty years-old patient, who finds her teeth too "yellow". We chose here a technique which apparently generates no sensitivities, especially for young patients. The White+™ technology includes a LED lamp allowing the activation of the bleaching gel's properties on the tooth surface. The patient's mouth preparation is less demanding: a good gingival isolation is enough as a condition to complete the treatment.

Picture 4b: Results at the end of the treatment (complete treatment duration at the dentist's chair = 45min)
We get good results in comparison with the starting case.

Features of the White+™ professional teeth whitening system

Conf. Univ. Dr. Anca Silvia Valceanu
Discipline of Dento-Facial Aesthetics, Faculty of Dental Medicine,
UMF „Victor Babeș” Timișoara

Introduction

Tooth whitening represents an important phase out of all the aesthetic rehabilitation treatments, which is a trend that's gaining more and more ground among today's adolescents and young adults.

All the hydrogen peroxide whitening products out on the market present basically the same working mechanism: an exothermal oxidation reaction which is mostly followed by teeth dehydration and post-op sensitivity.

The White+™ -professional system presents an original patent based on ecological active agents which maintain a neutral pH and optimal teeth hydration level; more than that, the LED lamp doesn't produce heat and presents 3 different operating programs for a better optimality of treatment results.

Purpose of study

The purpose of this study was to establish the efficiency of Meodental Company's White+™ - Professional System in depigmentation of dyschromia and to verify the presence or absence of post-op sensitivity in young patients.

Material and method

The whitening treatment targeted young teeth whose natural color was more yellowish. 8 arches have been treated in patients aged between 23-35 years old; all patients wanted a diminishment of dental chroma up to A1-B1. Every patient has been applied with White+ professional whitening gel in optimal isolation conditions; the gel was activated using White+ LED Lamp set on the corresponding program, according to the producers indications.

Results

The whitening treatments results was amazing after just one appointment whose duration didn't exceed 30 minutes. Dental sensitivity was spotted in just one case, in one fissured central incisive. Under no circumstances was a home treatment addition necessary.

Conclusions

Once familiarised with Meodental Company's White+™ - professional whitening system, it is easy to affirm that this presents a series of advantages: easy to use, occupying a small space in your practice, efficient with just one appointment, painless procedure which does not end up with tooth dehydration.

Key words: White+ professional, dental cosmetics, tooth whitening

1. Smile before treatment
2. Initial teeth color (Vitapan Classical)
- 3.- 6. The treatment sequence.
The treatment was performed with White + professional gel and with the White+ lamp progressive program for 30 min.
7. The color of the superior teeth after one treatment session
8. The bleaching treatment at the mandibular level had to be interrupted after 10 minutes due to sensitivity generated by the fissures existing at this level.
9. Whitening of the mandibular teeth continued at home using White+ Home gel during one week.
The image illustrates the tray applied at the inferior level.
10. Post-treatment smile
11. The satisfied patient for a successful result.

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To MEODENTAL
Dror Kanion

Dear Sirs,

Following a one year experience with Meodental White + professional, we find the kit most beneficial in several aspects:

1. Easy and time saving patient preparation – convenient both to the dentist and to the patient.
2. The bleaching result is visible after a 30 min. cycle of light exposure. In cases of extreme staining, adding a 15 minutes treatment of the "intensive" program improves results.
3. Non-painful procedure – very few patients reports of sensitivity, both during treatment and post-operative.
4. Long lasting results – combined with home bleach materials, the result is very stable.

The results are excellent and most agreeable to our patients.

Sincerely,

Dr. Galit Talmor

Clinic of Nitzan Bichacho D.M.D.

Bichacho Clinic - 14 Weizmann St. 64239, Tel-Aviv, ISRAEL
Tel: +972 3 6050060 Fax: +972 3 6024862 Mobile: +972 54 4534949
www.leadingdentalcenters.com

Professor and Head

R. E. Goldstein Center for Aesthetic Dentistry and Clinical Research
Prosthodontics Dep. Faculty of Dentistry, Hebrew University - Hadassah, Jerusalem

רוז ויצמן 14, תל אביב 64239 טל. 03-6054370 פקס. 03-6024862
14 Weizmann St., Tel Aviv 64239 Israel, Tel. 972-3-6050060, Fax. 972-3-6024862
E-Mail: bichacho-office@bichacho.net

Meodental at IDS 2013

meodental

meodental

www.meodental.com

meodental

Contact

For additional information about our products, please feel free to contact us:

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